

SAY SKIN  
AURORA<sup>+</sup>

Simple Manual

SAY SKIN  
AURORA<sup>+</sup>

**SAY SKIN** is  
the premium beauty care device brand  
that is designed to help achieve flawless skin.



SAY SKIN AURORA<sup>+</sup>

SAY SKIN A/S : 1661-6987

Copyright ©2020 Withnix Co., Ltd. All Rights Reserved

**MADE IN KOREA**

Stay beautiful forever with SAY SKIN

## # HOW TO USE THE DEVICE



Press the On/Off button long enough

Switch it On/Off

- Level 3
- Level 2
- Level 1



Press the Level button shortly. Each time you press the button, the level will go up. If you press the button at Level 3, it will go back to Level 1.

Choose Level

## # HOW TO CARE

STEP 01.




Clean the skin after washing your skin.

STEP 02.




Apply SAY SKIN super moisture cream.

STEP 03.

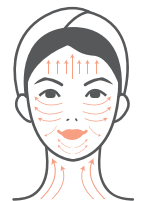

Care your skin from bottom to top as the picture.

## # HOW TO CHARGE THE DEVICE

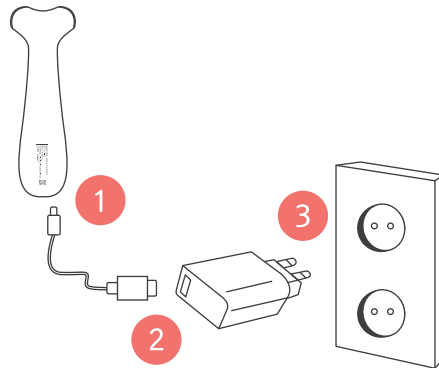

- 1) Connect the USB cable (Type C) to the main body of the device (**Fig. 1**)
- 2) Connect the USB cable (Type C) to the charging adapter (**Fig. 2**)
- 3) Put the adapter plug into the outlet (**Fig. 3**)

## # SAFETY PRECAUTIONS

- If you have any of the following symptoms or conditions, please consult your doctor before use.
    - 1) Rashes or excessive irritation appear on the skin.
    - 2) You need to take precautions because you are pregnant, in your menstrual period or in outpatient treatment.
    - 3) You have allergic diseases with highly sensitive skin.
    - 4) You have a serious infection caused by festering pimples.
  - People who use the following medical devices cannot use this product.
    - 1) Users of implantable electronic medical devices such as a pacemaker
    - 2) Users of electronic life support medical devices such as a heart lung machine.
    - 3) Users of electronic medical devices that are attached to the body such as ECG devices
  - If a person with sensitive skin uses this product, this may cause minor skin contact irritation. Stop use for two to three days if these symptoms appear, and then start using it again.
- ※ This is a skin care device, not a medical device, and does not have any medical efficacy or effect other than the beauty care effects.